




Waterbath

WNB 29

The combination of corrosion-resistant stainless steel, precise electronics and multiple temperature protection guarantees the highest level of safety in the laboratory.



On this page, you can find all the essential technical data on the Memmert water bath. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at myAtmoSAFE@memmert.com.

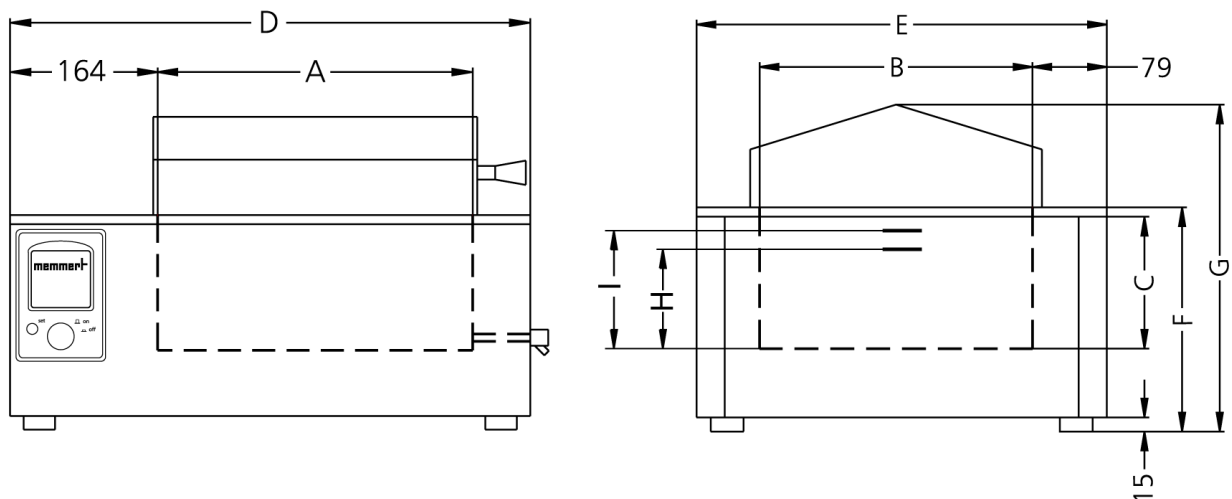

Temperature

Temperature range in °C	min. 5°C above ambient up to +95°C with additional boiling mode (+100°C)
-------------------------	--

Control of standard components

Temperature	1 Pt100 sensor class A in 4-wire-circuit
Timer	integrated digital timer from 1 min. to 99,59 hours for: ON continuous operation WAIT (delayed on for continuous and limited timed operation) HOLD

Safety

Temperature control	mechanical temperature limiter TB protection class 1 switching the heating off at approx. 30°C above max. temperature of the bath
Temperature control	in case of overtemperature due to failure, the heating is switched off at approx. 10°C above the set temperature (fixed value)
Autodiagnostic system	microprocessor PID-temperature controller with integrated autodiagnostic system with fault indicator

Heating concept

Heating Baths	corrosion-proof large-area heating on three sides
---------------	---

Standard equipment

Controller	digital display (LED) of set and actual temperature (0,1°C resolution) and of (remaining) programme time
------------	--

Stainless steel interior

Interior	easy-to-clean interior, made of stainless steel, reinforced by deep drawn ribbing, material no. 1.4301 (ASTM 304), laser-welded
Volume	29 l
Dimensions W x H x D in mm	$w_{(A)} \times h_{(B)} \times d_{(C)}$: 590 x 350 x 140 mm
Liquid level min.	(H) 105 mm
Liquid level max.	(I) 120 mm

Textured stainless steel casing

Dimensions	$w_{(D)} \times h_{(E)} \times d_{(F)}$: 818 x 516 x 238 mm
------------	--

Electrical data

Voltage 230 V, 50/60 Hz

Electrical load approx. 2400 W (during heating)

Packing/shipping data

the appliances must be transported upright

Customs tariff number 8419 8998

Country of origin Federal Republic of Germany

WEEE-Reg.-No. DE 66812464

Dimensions approx incl. carton B x H x T: 910 x 610 x 400 mm

Net weight approx. 24 kg

Gross weight carton approx. 31 kg

Standard units are safety-approved and bear the test marks

